


VIPER VACUUM STABILIZER

Soft and Secure Beating Heart Stabilization Technology


The Viper Vacuum Stabilizer™ combines Soft-Touch™ footpads for minimal tissue trauma and excellent stability for access to all coronary arteries.



**Soft-Touch
Rotating Footpads**

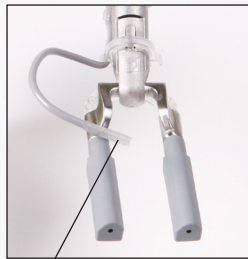

**Optimal Integrated
Blower-Mister**

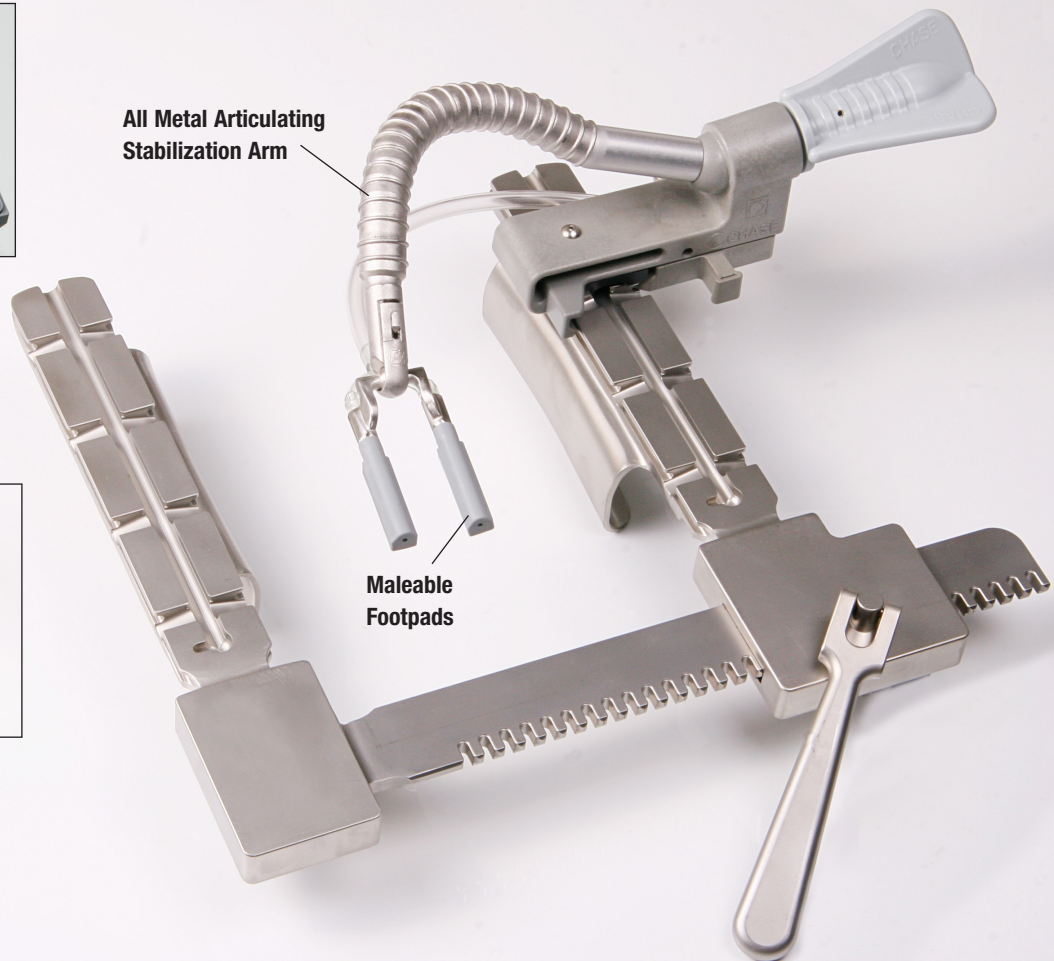

**All Metal Articulating
Stabilization Arm**

**Maleable
Footpads**

Viper Features

- *Soft-Touch footpads minimize tissue trauma*
- *Malleable footpads provide optimal positioning*
- *Rotating footpads thin epicardial fat*
- *Integrated Blower-Mister optional feature*
- *Stabilizer fits other sternal retractors*


VIPER VACUUM STABILIZER SYSTEM

The single-use Viper Vacuum Stabilizer with its unique footpad design and flexible material delivers superior stability and performance during beating heart procedures.



Soft-Touch Vacuum Footpad

- Raised-ridge cavities lock in vacuum
- Maintains suction on uneven cardiac contours

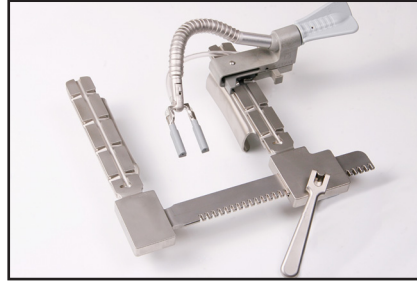

Retractor System

- Low profile rail system with disposable suture stays
- 3-sided mounting



Attached Blower-Mister

- Eliminates need for operator assistance
- "Hands free" operation

ORDERING INFORMATION

CODE	DESCRIPTION	UNITS PER CARTON
STB-50CVS	Viper Heart Stabilizer, C Mount	1
STB-5ACVS	Viper Heart Stabilizer, C Mount, Blower-Mister	1
STB-50MVS	Viper Heart Stabilizer, M Mount*	1
STB-5AMVS	Viper Heart Stabilizer, M Mount, Blower-Mister*	1
STB-50QVS	Viper Heart Stabilizer, Q Mount*	1
STB-5AQVS	Viper Heart Stabilizer, Q Mount, Blower-Mister*	1
FTV-1010	Single Foot Pedal for Stabilizer, Blower-Mister	1
VTS-0001S	Viper Tubing Set	1
RTC-1000	Sternal Retractor – Reusable	1
RTC-1020S	Suture Stays – Disposable	5 pairs

*We offer two additional base designs to mount on other sternal retractors.

TO PLACE ORDERS

Chase Medical Customer Service
 1-800-783-0378 or 972-783-7005
 Fax: 972-235-3446
 Email: cust_svc@chasemedical.com